



**QUEEN'S
UNIVERSITY
BELFAST**

Development and validation of an ultrasensitive fluorescence planar waveguide biosensor for the detection of paralytic shellfish toxins in marine algae

Meneely, J. P., Campbell, K., Greef, C., Lochhead, M. J., & Elliott, C. T. (2013). Development and validation of an ultrasensitive fluorescence planar waveguide biosensor for the detection of paralytic shellfish toxins in marine algae. *Biosensors and Bioelectronics*, 41(null), 691–697. <https://doi.org/10.1016/j.bios.2012.09.043>

Published in:

Biosensors and Bioelectronics

Document Version:

Peer reviewed version

Queen's University Belfast - Research Portal:

[Link to publication record in Queen's University Belfast Research Portal](#)

Publisher rights

This is the author's version of a work that was accepted for publication in *Biosensors and Bioelectronics*. Changes resulting from the publishing process, such as peer review, editing, corrections, structural formatting, and other quality control mechanisms may not be reflected in this document. Changes may have been made to this work since it was submitted for publication. A definitive version was subsequently published in *Biosensors and Bioelectronics*, VOL 41, 03/2013

General rights

Copyright for the publications made accessible via the Queen's University Belfast Research Portal is retained by the author(s) and / or other copyright owners and it is a condition of accessing these publications that users recognise and abide by the legal requirements associated with these rights.

Take down policy

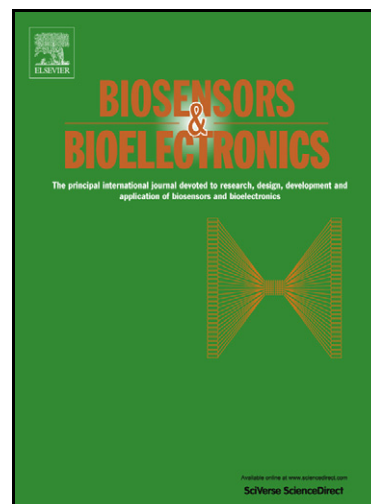
The Research Portal is Queen's institutional repository that provides access to Queen's research output. Every effort has been made to ensure that content in the Research Portal does not infringe any person's rights, or applicable UK laws. If you discover content in the Research Portal that you believe breaches copyright or violates any law, please contact openaccess@qub.ac.uk.

Open Access

This research has been made openly available by Queen's academics and its Open Research team. We would love to hear how access to this research benefits you. – Share your feedback with us: <http://go.qub.ac.uk/oa-feedback>

Development and validation of an ultrasensitive fluorescence planar waveguide biosensor for the detection of paralytic shellfish toxins in marine algae

Julie P. Meneely, Katrina Campbell, Charles Greef, Michael J Lochhead, Christopher T Elliott



www.elsevier.com/locate/bios

PII: S0956-5663(12)00649-5
DOI: <http://dx.doi.org/10.1016/j.bios.2012.09.043>
Reference: BIOS5439

To appear in: *Biosensors and Bioelectronics*

Received date: 31 July 2012
Revised date: 12 September 2012
Accepted date: 22 September 2012

Cite this article as: Julie P. Meneely, Katrina Campbell, Charles Greef, Michael J Lochhead and Christopher T Elliott, Development and validation of an ultrasensitive fluorescence planar waveguide biosensor for the detection of paralytic shellfish toxins in marine algae, *Biosensors and Bioelectronics*, <http://dx.doi.org/10.1016/j.bios.2012.09.043>

This is a PDF file of an unedited manuscript that has been accepted for publication. As a service to our customers we are providing this early version of the manuscript. The manuscript will undergo copyediting, typesetting, and review of the resulting galley proof before it is published in its final citable form. Please note that during the production process errors may be discovered which could affect the content, and all legal disclaimers that apply to the journal pertain.

Development and validation of an ultrasensitive fluorescence planar waveguide biosensor for the detection of paralytic shellfish toxins in marine algae

Julie P. Meneely^{a*}, Katrina Campbell^a, Charles Greef^b, Michael J Lochhead^b, Christopher T Elliott^a

^a*Institute of Agri-food and Land Use, School of Biological Sciences, Queen's University, Malone Road, Belfast, BT9 5BN, United Kingdom*

^b*MBio Diagnostics Inc, Boulder, Colorado, 80301, USA*

*Corresponding author at: Institute of Agri-food and Land Use, School of Biological Sciences, Queen's University, Malone Road, Belfast, BT9 5BN, United Kingdom

Tel: +44 (0) 2890 976562

Fax: +44 (0) 2890 976513

Email address: j.p.meneely@qub.ac.uk

ABSTRACT

Marine dinoflagellates of the genera *Alexandrium* are well known producers of the potent neurotoxic paralytic shellfish toxins that can enter the food web and ultimately present a serious risk to public health in addition to causing huge economic losses. Direct coastal monitoring of *Alexandrium* spp. can provide early warning of potential shellfish contamination and risks to consumers and so a rapid, sensitive, portable and easy-to-use assay has been developed for this purpose using an innovative planar waveguide device.

The disposable planar waveguide is comprised of a transparent substrate onto which an array of toxin-protein conjugates is deposited, assembled in a cartridge allowing the introduction of sample, and detection reagents. The competitive assay format uses a high affinity antibody to paralytic shellfish

toxins with a detection signal generated via a fluorescently labelled secondary antibody. The waveguide cartridge is analysed by a simple reader device and results are displayed on a laptop computer. Assay speed has been optimised to enable measurement within 15 min. A rapid, portable sample preparation technique was developed for *Alexandrium* spp. in seawater to ensure analysis was completed within a short period of time. The assay was validated and the LOD and CC β were determined as 12pg/mL and 20pg/mL respectively with an intra-assay CV of 11.3% at the CC β and an average recovery of 106%. The highly innovative assay was proven to accurately detect toxin presence in algae sampled from the US and European waters at an unprecedented cell density of 10cells/L.

Keywords: Paralytic shellfish toxins, fluorescence immunoassay, planar waveguide, *Alexandrium* spp, early warning, ultra sensitive.

1. Introduction

Phytoplankton is intrinsic to the marine ecosystem as they produce a significant proportion of the photosynthetic biomass that zooplankton and larger life forms rely on. Under specific conditions of high nutrient levels from agricultural, domestic and industrial discharges, these microalgae proliferate and produce harmful algal blooms (HABs, also known as “red tides”) (James et al., 2010). Various HAB species are toxic and can lead to poisoning of humans, fish and other wildlife via food web transfer (Batoréu et al., 2005). While many dinoflagellate genera cause HABs, one of the most important is the *Alexandrium* genus because of its global distribution, diversity and severity of the toxins produced (Churro et al., 2010). This genus is responsible for the production of saxitoxins (STX) that cause paralytic shellfish poisons (PSPs) throughout the world (Anderson et al., 2012).

PSP in humans is caused by the consumption of contaminated bivalves and other shellfish that have accumulated STX during filter feeding. The STX are potent neurotoxins that bind to the sodium channels

in nerve cell membranes disrupting normal depolarization. To date 57 STX analogues have been identified and are collectively referred to as paralytic shellfish toxins (PSTs) (Wiese et al., 2010), all exhibiting different toxicities.

Coastal monitoring for the presence and identification of marine toxic algae and the concentrations of PSTs produced is a significant challenge (Anderson, 2009). That said while monitoring is clearly required, the costs must be controlled. As a result, field deployable, sensitive, rapid qualitative tests will help provide protection of these areas.

Many methods of analysis exist for the measurement of the PSTs but generally these assays are designed for shellfish. In the absence of any regulations for PST concentrations in coastal waters, these techniques may not have the necessary sensitivity to provide early warning of these toxins. The reference method of analysis for shellfish is the mouse bioassay (MBA) (AOAC Official Method 959.08, 2005) however it is not sensitive enough as an early warning detection of PSTs in coastal waters. In addition to its technical shortcomings for water analysis, ethical considerations are becoming more important globally (Campbell et al., 2011a). The AOAC HPLC method (AOAC Official Method 2005.06., 2005) has found acceptance for screening, however sample preparation is lengthy and the interpretation of the results extremely complex requiring highly trained personnel (van den Top et al., 2011). Another HPLC method frequently used is that developed by Oshima (1995) however three different analyses are required to identify the PST derivatives. In recent times LC-MS/MS has emerged as a powerful tool for the analysis of PSTs (Jaime et al., 2001 and Dell'Aversano et al., 2005). While these methods allow for the accurate determination of these toxins in phytoplankton, they are laboratory based requiring skilled technicians and expensive equipment. More rapid immunological based methods have emerged over the past decade. One such example is biosensor assays based on surface plasmon resonance (Fonfria et al., 2007; Rawn et al., 2009; Campbell et al., 2010 and Haughey et

al., 2011), however again these have been designed for determination of the toxins in shellfish, are laboratory based and require a large capital outlay.

Many ELISAs have been applied to the detection of PSTs in shellfish samples but few have been tested for phytoplankton. They are specifically designed for screening, but are not truly portable requiring a reader. Often the sensitivity is high, but inaccurate results may occur as a result of the antibody specificity (Usleber et al., 2001; Dubois et al., 2010). One approach to alleviate the problem of underestimation of results is to combine two antibodies into the test format as described by Campbell et al (2011b) or employ two assays (Chu et al, 1996). Several commercial ELISA kits exist for the screening of PSTs, however only one has been applied to the determination of STX in water samples, the Abraxis Saxitoxin (PSP) ELISA kit. An alternative immunological format is that of Lateral Flow Devices (LFDs) which are designed to be used in field conditions and give a qualitative result, but again they have only been tested in shellfish.

A receptor binding assay (Doucette et al, 1997) and cytotoxicity/cell culture methods have been developed (Garthwaite, 2000) however while useful for research they may be very difficult to introduce for monitoring purposes. Field deployable techniques such as a neuronal network biosensor (Kulagina et al., 2006) and the use of synthetic resins to adsorb PSTs (Rodríguez et al, 2011) have been published. This study investigated the use of fluorescence immunoassay planar waveguide as an early warning detection system for PSTs produced by the *Alexandrium* spp. The assay has been shown to be rapid, easy-to-use and is ultra sensitive allowing detection of these toxins at very low algal cell densities. It can be applied in the field and thus raises the possibility of providing, for the first time, a fit-for-purpose monitoring tool for coastal regions often affected by HABs.

2. Materials and methods

2.1. Apparatus

A Minimix standard shaker was purchased from Merris Engineering, Maidenhead, UK. The SnapEsi[®] LS System and microarray cartridges were supplied by MBio Diagnostics Inc. (Boulder, USA). BagPage[®] + 80 (full surface micro-perforated filter blender bags, filter 280 µm) were supplied by Interscience, France.

2.2. Reagents

All chemicals were of Analar grade unless otherwise stated. Phosphate buffered saline (PBS) tablets, Tween 20, bovine serum albumin (BSA), dimethylsulphoxide, boric acid, tris (hydroxymethyl) aminomethane, 2,2'-(ethylenedioxy)bis(ethylamine) (jeffamine, MW 148.2) and HPLC grade water and acetonitrile were purchased from Sigma-Aldrich, United Kingdom.

Polyclonal antibodies against saxitoxin (QUB 5) and neosaxitoxin (QUB 7), (1mg/mL) and toxin protein conjugates, saxitoxin-ovalbumin and neosaxitoxin-ovalbumin (1mg/mL) were prepared in-house. Alexa Fluor[®] 647 goat anti-rabbit IgG, (2mg/mL) and Alexa Fluor[®] 647 carboxylic acid, succinimidyl ester, (1mg), were purchased from Invitrogen Ltd, Paisley, Scotland.

PST standards saxitoxin (STX), decarbamoyl saxitoxin (dcSTX), neosaxitoxin (NEO), decarbamoyl neosaxitoxin (dcNEO), gonyautoxins (GTX1/4 and GTX 2/3), decarbamoyl gonyautoxins (dcGTX 2/3), GTX5 (B1) and C1/2 standard solutions were obtained from the National Research Council of Canada, (NRC), Halifax, Canada. All standards were diluted in PBS (pH 7.2) to prepare working standards for the immunoassay and HPLC analysis. For the preparation of the toxin protein conjugates, STX and NEO were also purchased from the NRC.

Strains of *Alexandrium tamarense* (Lebour) Balech, CCAP 1119/20, 1119/27 and 1119/32 were purchased from the Culture Collection of Algae and Protozoa, Scottish Marine Institute, Oban, Scotland (CCAP). Two further strains, *Alexandrium fundyense* (Lebour) Balech, F14 and GTCA28 were supplied by the School of Marine Sciences, University of Maine, USA. L1 medium for marine dinoflagellates was

purchased from the CCAP. Fresh natural seawater was obtained from Strangford Lough, Co. Down, Northern Ireland.

2.3. The assay cartridge and instrument

Assays were performed using the SnapEsi[®] LS System from MBio Diagnostics Inc. The SnapEsi[®] System is a multiplexed, fluorescence immunoassay platform in development for clinical diagnostics and life sciences applications. It combines single-use disposable cartridges with a simple reader, and is described in detail elsewhere (Lochhead et al., 2011). The system utilizes a novel implementation of planar waveguide technology. Planar waveguides for immunoassays have been the subject of several technical reviews (Herron et al., 2005; Ligler., 2009; Mukundan et al., 2009). Briefly, a light source (typically a laser) is directed into a waveguide substrate where it propagates by total internal reflection (TIR) at the interface between the high index of refraction waveguide (glass or plastic) and the surrounding medium (air or aqueous solution). At the TIR interface, an evanescent field is generated that decays exponentially into the aqueous medium. This provides localized illumination of the assay surface, minimising effects of the bulk solution, line-of-sight, light scattering, etc.

The key innovation of the MBio Diagnostics planar waveguide system is an optical design that provides an easy-to-use, reproducible, all plastic disposable cartridge that avoids the cost and alignment issues that have limited adoption of planar waveguide technology in high volume applications.

2.4. Microarray construction

Toxin-protein conjugates were printed to chemically activated waveguides (Grainger et al., 2007) at defined concentrations as depicted in Figure 1. Fluorescently labelled protein conjugate was printed as a spatial orientation marker, and anti-species antibodies served as procedural controls. Arrays were printed with a Bio-Dot AD3200 robotic arrayer equipped with a Bio-Jet print head dispensing 20nL

droplets. Spot diameters were approximately 0.5mm and the arrays were printed on a grid with 1mm centres. For the final array, each toxin-protein conjugate was printed in replicates of four. After printing, the waveguide arrays were rinsed with a protein-based blocking agent and dried by centrifugation.

Printed waveguides were assembled into an injection moulded cartridge using die cut double sided adhesive, forming a 5mm wide fluidic channel with a volume of approximately 30 μ L. The cartridge inlet port provides a reservoir for introduction of assay fluids. All assay fluids stay on-board the cartridge (minimizing biohazard), with total waste volume up to 1 mL irreversibly trapped in an absorbent pad. The cartridge has been engineered to provide reproducible, passive fluid flow, driven by a combination of capillary action and hydrostatic pressure.

2.5. Toxin protein conjugates

STX was conjugated to the carrier protein, keyhole limpet haemocyanin (KLH) via a modification of the Mannich reaction incorporating jeffamine as a spacer compound and polyclonal antibodies (PAb) were prepared in a similar manner to that reported previously for bovine serum albumin (Campbell et al., 2007). The antibody was purified by ammonium sulphate precipitation followed by dialysis in 0.15M saline and the protein concentration was determined by A_{280} to be 8.3mg/mL.

For the microarray, neosaxitoxin was conjugated to the carrier protein ovalbumin using the same conjugation chemistry as above in the following ratios but omitting the jeffamine linker. Ovalbumin (1mg), neosaxitoxin (146 μ g) and formaldehyde (4.5 μ L) were reacted followed by dialysis over 24 h in 3 x 4 L of 0.15M saline solution.

2.6. Labelling the primary antibody with Alexa Fluor® 647 carboxylic acid, succinimidyl ester

In brief, 50µg of the saxitoxin antibody was placed in an Amicon UltraCel 30K membrane centrifugal filter and 100µL 50mM borate buffer, pH 8.0 added. This was centrifuged twice at 13,400 x g for 10 min, the protein solution removed and Alexa Fluor® 647 carboxylic acid, succinimidyl ester added (20-30 equivalents of dye to protein). These were mixed vigorously for 2 h at room temperature, shielded from direct sunlight. The reaction was quenched by the addition of an equal volume of 0.1M tris buffer, pH 8.5 and the contents transferred to another centrifugal filter. The volume was made up to 400µL with PBS containing 0.05% Tween 20 and centrifuged at 13,400 x g for 10 min. The supernatant was discarded and 400µL of the above buffer added and centrifuged as before. This process was repeated 3 times until the supernatant was colourless. The labelled antibody was stored at -20°C until use.

2.7. Antibody specificity

Antibody cross-reactivity was determined in assay buffer (PBS, pH 7.2, 1% BSA, 0.05% w/v Tween 20) spiked with structurally related analytes. The calibration points were 0, 0.02, 0.04, 0.08 and 1.6 ng/mL. The fluorescence signal obtained for each calibration point ($n = 4$) was normalized by calculating the percentage of bound antibody for each, relative to the zero calibration signal observed. The mean values were plotted using the non-linear four parameter equation using BIA evaluation software. The midpoint concentration (IC_{50}) of each toxin curve was determined and the cross-reactivity calculated as a percentage relative to STX.

2.8. Immunoassay procedure

The microarray cartridge was pre-conditioned with assay buffer (150µL) for 10 min. Antibody (100µL) and STX standard/sample (100µL) were mixed together and 150µL applied to the cartridge for 15 min

followed by two sample rinses using assay buffer (2 x 150µL). The labelled antibody (Alexa Fluor 647 goat anti-rabbit IgG at a concentration of 10µg/mL; in assay buffer) was applied to the cartridge and incubated for 10 min. Two label rinse steps were performed as before and the data recorded on the SnapEsi® System. This assay procedure required 45 min to complete. Optimal working dilutions of antibodies and toxin-protein conjugates were determined experimentally by the checkerboard titration method using this protocol.

2.9. Algal culturing

L1 medium was prepared according to the manufacturer's guidelines in filtered seawater.

All glassware and media were sterilized prior to use and subsequent protocols performed in a UV3 HEPA PCR cabinet. The cultures were grown at 20 °C ±

(5 °C) with a 12:12 hour light: dark cycle with a photon flux density of approximately 170 µmol/m²/s in L1 medium. Batch cultures of each strain were inoculated from a stock culture that had reached the late exponential/early stationary phase and the cultures maintained in duplicate under the conditions described above.

2.10. Algal cell lysis and extraction of PSTs

Chemical and mechanical methods were investigated. Briefly, 10/15mL of culture was filtered using Millipore hydrophilic mixed cellulose esters filters (0.45µm). The algal cells were washed off the filter prior to reconstitution in 0.9/1.5mL of PBS. For chemical lysis 100µL copper sulphate/Lysozyme was added to the algal cells for 10 min. Mechanical lysis employed a metal seam roller to crush the cells in sample bags for 1.5 min. Concentrated cultures were also placed into tubes and whisked for 2 min with a hand blender. Cell counts were performed prior to and after each treatment as follows. Triplicate 1mL volumes of each culture were collected and fixed with Lugol's iodine solution. 1mL was transferred

to a Sedgewick-Rafter counting chamber for enumeration. Cell counts were performed manually at x 100 magnification on an inverted light microscope.

The 'gold' standard mechanical lysis method used was based on a modified version of glass bead beating described by Devlin et al (2011). Briefly 10mL of culture was concentrated by centrifugation, the supernatant discarded and the cells re-suspended in PBS. Glass beads were added and the samples placed in a Minimix standard shaker for 10 min. Cell enumeration was performed prior to extraction to enable the determination of total PST per cell in addition to the culture toxin profiles. Extracts from this method and the most efficient portable method were analysed by HPLC.

2.11. HPLC analysis of algal extracts

All extracts were analyzed to characterize the toxin profiles and quantify the toxin concentrations in each. The PSP toxin concentrations from the standard and portable lysis method selected were compared. The official AOAC method was used (AOAC Official Method 2005.06., 2005).

2.12. Assay performance

Limit of detection (LOD) for the assay was determined by the analysis of a non-toxic strain of *Alexandrium tamarense* (Lebour) Balech, CCAP 1119/20 ($n = 20$) and calculated as the mean of the measured fluorescence signal minus 3 times the standard deviation of the mean (mean - 3sd). The detection capability CC β was evaluated by fortification of the non-toxic strain at the concentration 20pg/mL and analyzed according to the assay protocol. The final assay calibration curve was 0, 0.02, 0.04 and 0.08ng/mL prepared by spiking an extract of the non-toxin producing strain at a cell density of 50cells/L. To determine the sensitivity of the assay with respect to cell numbers, both non-toxic and

toxic strains of *Alexandrium* were subjected to analysis at 50, 10 and 5 cells/L ($n = 6$ for each strain). Cell counts were performed for each followed by dilutions to the required densities prior to analysis.

3. Results and discussion

The competitive assay format uses a high affinity antibody to paralytic shellfish toxins with a detection signal generated via a fluorescently labelled secondary antibody.

The developed assay was designed to provide early warning of PSTs in coastal waters and provides rapid, qualitative results. As such, the sensitivity, specificity, LOD and CC β were evaluated to ensure the method was fit-for-purpose. The reproducibility of the printed waveguides was ascertained by performing both buffer and matrix curves over a number of days and using different print batches (see the supporting information section SI1). The standard deviations associated with these were all less than 5 with CV values ranging from 8% to 68% (25 and 250 pg/mL) for buffer (2 batches) (Table SI1). For the matrix curves (3 batches) the standard deviations were less than 6 and the CV values ranged from 0.6% to 56% (20 and 80 pg/mL) (Table SI2). The IC₅₀ values for the calibration curves were consistent at 24 pg/mL. The buffer and matrix curve overlays have been shown in the supplementary section (Figures SI1 and SI2).

The checkerboard titration method highlighted the QUB 5 polyclonal antibody (at a dilution of 1/8,000) provided the most sensitive results with the NEO-ovalbumin conjugate (50 μ g/mL). Under buffer conditions an IC₅₀, the concentration of toxin required to cause 50% inhibition of the antibody binding, of 24 pg/mL was demonstrated.

One of the main aims during the assay development was to achieve a short assay time. Initially the assay procedure was 45 min and considered too lengthy for a rapid test. To accomplish a faster analysis time, all washes were removed and the incubation times reduced to 7.5 min each for the

antibody/sample followed by the detection antibody. The IC_{50} was determined as 25pg/mL indicating no adverse effects yet a marked improvement in assay time was gained.

Further studies were performed to ascertain if a one-step protocol was possible. Two approaches were investigated, the first mixing the antibody, sample and secondary antibody together and incubating for 15 min prior to reading the data and the second was labelling the primary antibody with the fluorescent tag and applying the sample/antibody mixture for 15 min prior to reading. The first approach failed as very little binding was exhibited and the second, while it proved more successful showed a significant drop in sensitivity with the IC_{50} being calculated as 190pg/mL.

The antibody specificity was evaluated in assay buffer and the cross-reactivities were calculated as 100%, 25%, 28%, 8%, 49%, <0.1%, 7%, 33% and 4 % for STX, C1/2, dcGTX2/3, dcNEO, dcSTX, GTX1/4, GTX2/3, GTX5 and NEO respectively. The antibody was very specific for STX. Some cross-reactivity was exhibited towards dcSTX with less specificity towards dcGTX2,3 and dcNEO but this was not of particular relevance since decarbamoyl toxins are not commonly found in toxic dinoflagellates (Etheridge., 2010). Of significance was the low cross-reactivity of the antibody with the more potent STX analogues NEO and GTX1/4 and if present in the sample then ultimately the test could underestimate the results. These results are fairly consistent with other immunoassay publications (Usleber et al., 2001 and Dubois et al., 2010) however by use of a “cocktail” of antibodies, as described by Campbell et al (2011b) or two immunoassays (Chu et al., 1996), these difficulties may be overcome. Alternatively, the multiplexed nature of the microarray format could allow multiple separate tests to be run, enabling a comprehensive toxin profile assay.

An innovative approach was taken to develop a simple, rapid sample preparation technique compatible with the immunoassay. Importantly one of the necessary requirements of this assay was that it would be field deployable. The study was conducted using 3 cultures of *Alexandrium tamarense* including F14, GTCA28 and CCAP 1119/20. Chemical and mechanical lysis techniques were trialled as described and

the most efficient method chosen was initially based on cell counts and percentage lysis rates. Careful selection of the extraction buffer was fundamental, not only to ensure efficient extraction of the toxins but also to avoid any interference with the immunological processes. Instead of using acetic acid or hydrochloric acid as described by Ravn et al (1995), the authors opted for PBS, pH 7.2 (Devlin et al, 2011) to avoid adverse affects on the immunoassay. The results (Fig. 2) demonstrated the chemical lysis techniques investigated were not suitable with lysis rates of $\leq 54\%$. Physical disruption of the cell wall proved to be successful following crushing of the phytoplankton contained in the BagPage® + 80 sample bag. A lysis rate of $99.6 \pm 0.2\%$ was calculated for this method. Removal of the membrane from this system yielded an $83 \pm 7.4\%$ lysis rate while using a general laboratory sample bag or the hand blender provided lysis rates of only $54 \pm 2.8\%$ and $18 \pm 8.7\%$ respectively. Therefore the BagPage® + 80 and roller were selected. For analysis a filtration step was included after lysis to remove cell debris. The results observed using the new protocol were compared against the modified bead beating method (Devlin et al, 2011), by determining cell counts before and after lysis. Results using the modified bead beating method showed 100% lysis for all three cultures while $99 \pm 0.3\%$, $99 \pm 0.1\%$ and $98 \pm 0.2\%$ lysis was observed for cultures CCAP 1119/20, GTCA28 and F14 respectively using the BagPage® + 80 and roller. This indicated that the rapid portable method was highly satisfactory. Further investigations to compare the toxin concentrations using both the laboratory-based and portable methods were completed prior to validation of the method to confirm agreement. This was performed by repeating the tests as described and measuring the toxin concentrations by HPLC.

Table 1 highlights the results obtained from this trial. The total concentration of PSTs calculated using the bead beating method was taken as 100% and those calculated using the portable sample preparation method were measured against this as a reference. The recovery for the *Alexandrium* strains ranged from 71% to 91%, a highly acceptable result. In addition the total toxin concentrations

per cell were calculated for each algal strain, the highest level being produced by *Alexandrium fundyense* (F14) at the concentration 14.4 pg/cell.

The LOD of the assay was determined as described and a value of 12pg/mL was calculated. The lowest concentration of STX (CC β) that could be accurately detected was assessed by fortifying blank samples at 20pg/mL and the results presented in Fig. 3. The mean fluorescence signals for each sample were divided by the signal generated by the zero standard and multiplied by 100. The average normalized fluorescence signal $\pm$ 2sd were calculated for each population. No overlap was observed between them, the lower limit exhibited for the blank population was 80% while the upper limit for the spiked population was 68%. The Student's T Test was applied to the two sets of data and the p-value calculated as 3.9×10^{-19} , confirming the populations were significantly different. The intra-assay variation and extraction recovery were determined by the analysis of the samples spiked at 20pg/mL. The coefficient of variation (% CV) was calculated as 11.3% and the recovery rates ranged from 86% - 123%. Analysis of the *Alexandrium* strains at various cell densities allowed assessment of assay sensitivity. No guidelines exist for the measurement of toxin in algal cells however shellfish may begin to show signs of toxicity at about 100cells/L, and above 1000cells/L the toxin level in mussels may be above the regulatory limit (Blasco et al., 2003). Therefore, the ability to detect toxin in 50cells/L would allow significant early warning for potential contamination of shellfish. To assess if the assay could accurately determine if a sample was positive or negative at a cell density of 50cells/L, the cultures were diluted, cells lysed and the assay performed ($n = 6$) for each strain and the results are detailed in Table 2. In addition, the performance of the assay was investigated at lower cell densities of 10cells/L and 5cells/L. Total toxin concentrations determined by HPLC for the cell densities investigated were found to be greater than 20pg for cultures GTCA28, F14, CCAP 1119/27 and CCAP 1119/32. For CCAP 1119/20 (non-producer), the concentrations were less than 20pg. At 50cells/L and 10cells/L, the assay performance was exceptional, distinguishing between the negative and positive cultures. However at a

cell density of 5cells/L, although all the toxin producing algal strains were correctly identified, 2 samples of the non-toxin producing strain gave false positive results. Thus it can be concluded that the assay can perform very satisfactorily at detecting PST production in as low as 10cells/L and at this cell density no false negatives were observed during the validation study. Examining STX only, the assay performed well (Table 3). At 50cells/L there was complete agreement between the HPLC results and those obtained using the portable assay. At 10cells/L, the STX concentrations in the toxin-producing cultures by HPLC were all found to be below 20pg, while the portable assay found these to be positive, i.e. above the cut-off of 20pg. This is not surprising as the portable assay, although highly specific towards STX, will still pick up other toxins present due to the cross-reactivity profile of the antibody. With low levels of C1/2, GTX2/3, NEO and GTX5 in these algal samples, the portable assay will detect these and display positive results.

4. Conclusions

Monitoring of coastal waters for PSTs produced by marine dinoflagellates is of paramount importance to provide an early warning of risk. Rapid, on-site qualitative assays are advantageous in this situation because they are cost-effective and the generation of results is rapid. In this study a rapid, easy-to-use, portable fluorescent planar waveguide immunoassay has been developed and validated to meet these requirements. In conjunction with a novel sample preparation, the total time required to test a sample is approximately 20 minutes. The lowest detectable limit of the assay was in the low pg/mL range making it an ultra sensitive method. The assay has been designed to detect toxin levels in *Alexandrium* spp. known to have global distribution however the method is also likely to be applicable to other PST producing organisms in both the marine and freshwater environments.

Acknowledgements

The authors thank the Department of Employment and Learning via the All Island Research Programme for funding the Project: BEACONS: Biosafety for Environmental Contaminants using Novel Sensors (Reference, USI 001). We also thank Dr Laurie Connell and Amber Bratcher of the University of Maine for providing *Alexandrium* cultures and Dr Gregory Doucette, NOAA/National Ocean Service for his advice during the study.

References

- Anderson, D.M., 2009. Ocean Coast. Manage. 52, 342-347
- Anderson, D.M., Alpermann, T.J., Cembella, A.D., Collos, Y., Masseret, E., Montresor, M., 2012. Harmful Algae. 14, 10-35
- Batoréu, M.C.C., Dias, E., Pereira, P., Franca, S., 2005. Environ. Toxicol. Pharmacol. 19, 401-406
- Blasco, D., Levasseur, M., Bonneau, E., Gelinas, R., Packard, T., 2003. Sci. Mar. 67, 261-278
- Campbell, K., Haughey, S.A., van den Top, H., van Egmond, H., Vilarino, N., Botana, L.M., Elliott, C.T., 2010. Anal. Chem. 82, 2977-2988
- Campbell, K., Rawn, D.F.K., Niedzwiadek, B., Elliott, C.T., 2011b. Food Addit. Contam. A. 28, 711-725
- Campbell, K., Stewart, L.D., Doucette, G.J., Fodey, T.L., Haughey, S.A., Vilarino, N., Kawatsu, K., Elliott, C.T., 2007. Anal. Chem. 79, 5906-5914
- Campbell, K., Vilariño, N., Botana, L.M., Elliott, C.T., 2011a. TrAC. 30, 239-253
- Chu, F., Hsu, K., Huang, X., Barrett, R., Allison, C., 1996. J. Agric. Food Chem. 44, 4043-4047
- Churro C., Alverca, E., Sam-Bento, F., Paulino, S., Figueira, V.C., Bento, A.J., Prabhakar, S., Lobo, A.M., Calado, A.J., Pereira, P., 2010. J. Appl. Phycol. 21, 429-442.
- Dell'Aversano, C., Hess, P., Quilliam, M.A., 2005. J. Chromatogr. A. 1081, 190-201

- Devlin, R.A., Campbell, K., Kawatsu, K., Elliott, C.T., 2011. *Harmful Algae*. 10, 542-548
- Doucette, G.J., Logan, M.M., Ramsdell, J.S., Van Dolah, F.M., 1997. *Toxicon*. 35, 625-636
- Dubois, M., Demoulin, L., Charlier, C., Singh, G., Godefroy, S.B., Campbell, K., Elliott, C.T., Delahaut, P., 2010. *Food Addit. Contam. A*. 27, 859-868
- Etheridge, S.M., 2010. *Toxicon*. 56, 108-122
- Fonfria, E.S., Vilarino, N., Campbell, K., Elliott, C., Haughey, S.A., Ben-Gigirey, B., Vieites, J.M., Kawatsu, K., Botana, L.M., 2007. *Anal. Chem*. 79, 6303-6311
- Garthwaite, I., 2000. *Trends Food Sci. Technol*. 11, 235-244
- Grainger, D.W., Greef, C.H., Gong, P., Lochhead, J.L., 2007. *Methods Mol Biol*. 381, 37-57
- Haughey, S.A., Campbell, K., Yakes, B.J., Prezioso, S.M., DeGrasse, S.L., Kawatsu, K., Elliott, C.T., 2011. *Talanta*. 85, 519-526
- Herron, J.N., Wang, H-K., Tan, L., Brown, S.Z., Terry, A.H., Durtschi, J.D., Simon, E.M., Astill, M.E., Smith, R.S., Christensen, D.A., 2006. Planar waveguide biosensors for point-of-care clinical and molecular diagnostics in Thompson R.B. (Ed.), *Fluorescence sensors and biosensors*. Taylor & Francis CRC Press, Boca Raton, FL, pp. 283-332.
- Humpage, A.R., Magalhaes, V.F., Froscio, S.M., 2010. *Anal. Bioanal. Chem*. 397, 1655-1671
- Jaime, E., Hummert, C., Hess, P., Luckas, B., 2001. *J. Chromatogr. A*. 929, 43-49
- James, K.J., Carey, B., O'Halloran, J., van Pelt, F.N.A.M., Skrabakova, Z., 2010. *Epidemiol. Infect.* 138, 927-940
- Kulagina, N., Mikulski, C., Gray, S., Ma, W., Doucette, G., Ramsdell, J., Pancrazio, J., 2006. *Environ. Sci. Technol*. 40, 578-583
- Ligler, F.S., 2009. *Anal. Chem*. 81, 519-526

- Lochhead, M.J., Todorof, K., Delaney, M., Ives, J.T., Greef, C., Moll, K., Rowley, K., Vogel, K., Myatt, C., Zhang, X., Logan, C., Benson, C., Reed, S., Schooley, R.T., 2011. *J. Clin. Microbiol.* 49, 3584-3590
- Mukundan, H., Anderson, A.S., Grace, W.K., Grace, K.M., Hartman, N., Martinez, J.S., Swanson, B.I., 2009. *Sensors*. 9, 5783-5809
- Official Methods of Analysis of the Association of Official Analytical Chemists., 2005. 18th ed.; AOAC INT. 88-96
- Oshima, Y., 1995. *J. AOAC Int.* 78, 528-532
- Ravn, H., Anthoni, U., Christophersen, C., Nielsen, P., Oshima, Y., 1995. *J. Appl. Phycol.* 7, 589-594
- Rawn, D.F.K., Niedzwiadek, B., Campbell, K., Higgins, H.C., Elliott, C.T., 2009. *J. Agric. Food Chem.* 57, 10022-10031
- Rodríguez, P., Alfonso, A., Turrell, E., Lacaze, J., Botana, L.M., 2011. *Harmful Algae*. 10, 447-455
- Usleber, E., Dietrich, R., Burk, C., Schneider, E., Martlbauer, E., 2001. *J. AOAC Int.* 84, 1649-1656
- van den Top, H.J., Elliott, C.T., Haughey, S.A., Viarino, N., van Egmond, H.P., Botana, L.M., Campbell, K., 2011. *Anal. Chem.* 83, 4206-4213
- Wiese, M., D'Agostino, P.M., Mihali, T.K., Moffitt, M.C., Neilan, B.A., 2010. *Mar. Drugs*. 8, 2185-2211

1119/20	beating									
	BagPage®									
	+ 80	0	0	0	0	0	0	0	100	0
CCAP	Bead									
1119/27	beating	100	796	113	86	6	223	1324	100	9.7
	BagPage®									
	+ 80	74	663	106	58	8	183	1092	82	8.2
CCAP	Bead									
1119/32	beating	0	782	26	21	0	92	921	100	9.5
	BagPage®									
	+ 80	0	607	33	19	0	106	765	83	7.9

Table 2. Qualitative assessment of the portable assay at *Alexandrium* spp. cell densities of 50, 10 and 5 cells/L.

<i>Alexandrium</i> Strain (n = 6)	50 Cells		10 Cells		5 Cells	
	Negative	Positive	Negative	Positive	Negative	Positive
<i>Alexandrium tamarense</i> 1119/20	6	0	6	0	4	2
<i>Alexandrium tamarense</i> 1119/27	0	6	0	6	0	6
<i>Alexandrium tamarense</i> 1119/32	0	6	0	6	0	6
<i>Alexandrium fundyense</i> GTCA28	0	6	0	6	0	6
<i>Alexandrium fundyense</i> F14	0	6	0	6	0	6

Negative refers to concentrations <20pg/cell, Positive refers to concentrations >20pg/cell

Table 3. Performance of the portable assay for the qualitative determination of STX in *Alexandrium* spp. at 50 and 10cells/L.

<i>Alexandrium</i> Strain	50 Cells		10 Cells	
	HPLC	Portable Assay	HPLC	Portable Assay
	(STX pg/50 cells)	Negative/Positive	(STX pg/10 cells)	Negative/Positive
Alexandrium tamarens 1119/20	0	Negative	0	Negative
Alexandrium tamarens 1119/27	82	Positive	16	Positive
Alexandrium tamarens 1119/32	48	Positive	10	Positive
Alexandrium fundyense GTCA28	76	Positive	15	Positive
Alexandrium fundyense F14	93	Positive	19	Positive

Negative refers to concentrations <20pg/cell, Positive refers to concentrations >20pg/cell

Supporting Information

Development and validation of an ultrasensitive fluorescence planar waveguide biosensor for the detection of paralytic shellfish toxins in marine algae

Julie P. Meneely^{a*}, Katrina Campbell^a, Charles Greef^b, Michael J Lochhead^b, Christopher T Elliott^a

*^aInstitute of Agri-food and Land Use, School of Biological Sciences, Queen's University, Malone Road,
Belfast, BT9 5BN, United Kingdom*

^bMBio Diagnostics Inc, Boulder, Colorado, 80301, USA

Accepted manuscript

SI Procedure: Assessment of the reproducibility of the printed waveguides.

The reproducibility of the printed waveguides was determined as follows. Under buffer conditions, two batches of waveguides were tested on four separate occasions over a period of several months ($n = 16$) for each calibration point whereas matrix calibration curves were repeated three times using three different batches over three days ($n = 12$). The fluorescence signal for each calibration solution was normalized by calculating the percentage of bound ligand (antibody) for each calibration point relative to the zero calibration fluorescence signal, the mean values plotted and the results indicate that the print procedure was both repeatable and reproducible for both buffer (Figure SI1) and matrix (Figure SI2). In addition, the results have been tabulated (buffer curves, Table SI1 and matrix curves Table SI2), highlighting the standard deviations and %CVs for each. For pg levels of toxin, the CV values are acceptable and the highest CV values occur at the high concentrations, i.e. where the calibration curves flatten out and almost 100% inhibition is observed. Please note the buffer calibration curves were performed during the early assay development stages and the points used were different from the final assay. For buffer the calibration points used were 0, 0.025, 0.05 and 0.25ng/ml while for matrix they were changed to 0, 0.02, 0.04 and 0.08ng/ml.

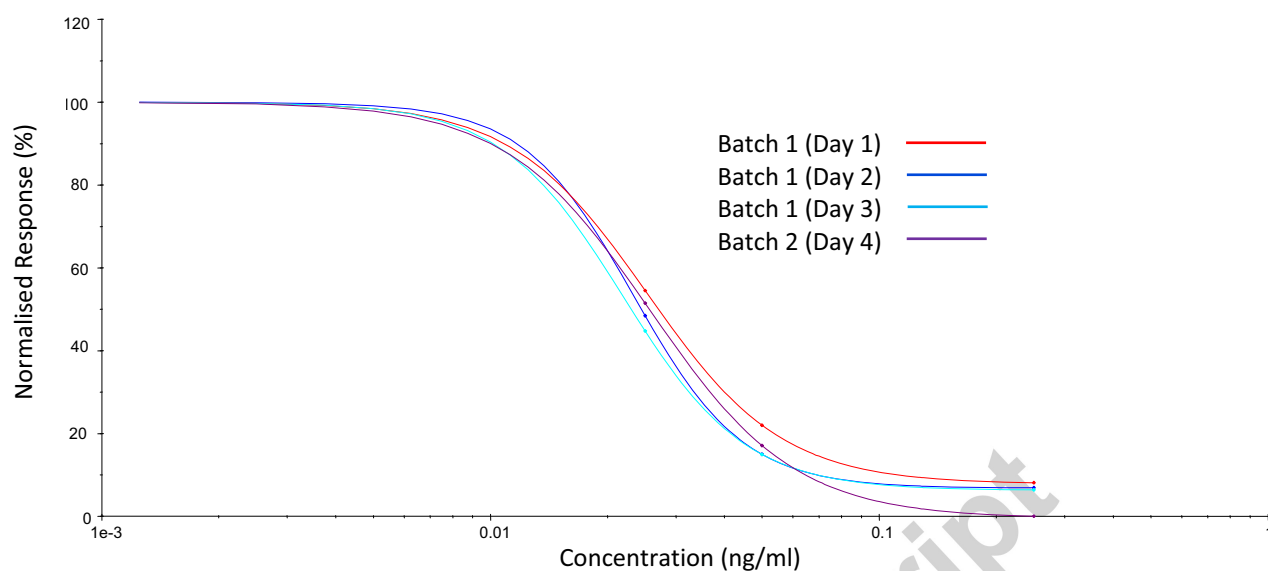


Fig SI1. Buffer calibration curves performed using 2 batches of printed microarrays over a period of several months (September to March).

*Colour reproduction on the Web.

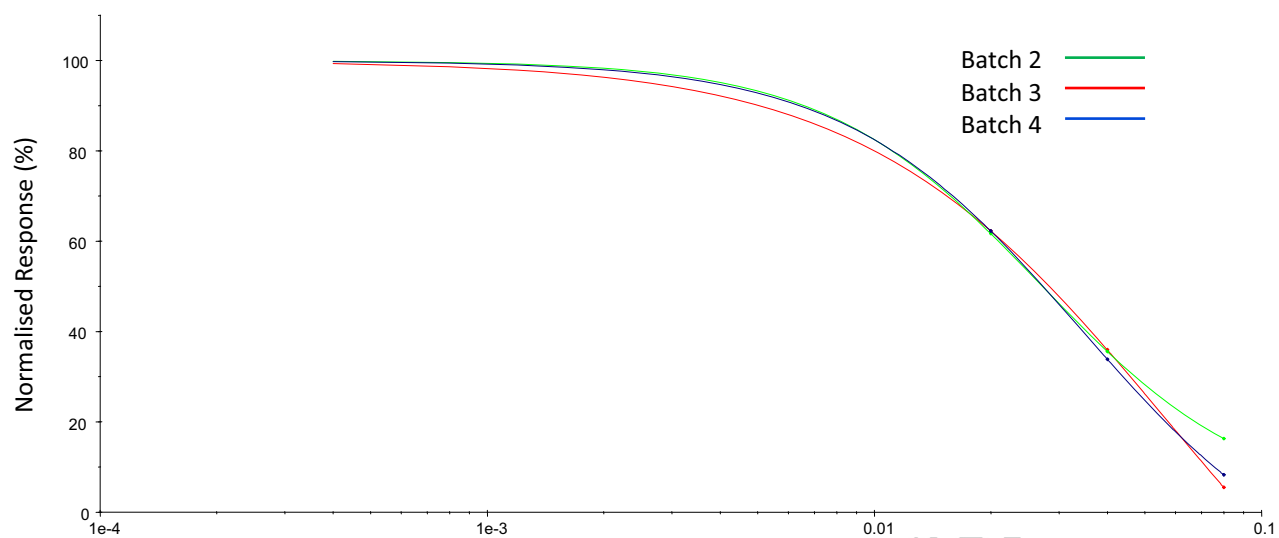


Fig SI2. Matrix calibration curves performed using 3 batches of printed microarrays over 3 different occasions.

*Colour reproduction on the Web.

Table SI1. Assessment of waveguide reproducibility by comparing buffer calibration curves on 2 batches (4 occasions) over several months.

Concentration	Day 1	Day 2	Day 3	Day 4	Mean	sd	CV
(ng/mL)	Normalised	Normalised	Normalised	Normalised			(%)
	Signal	Signal	Signal	Signal			
	(n = 4)	(n = 4)	(n = 4)	(n = 4)			
0	100	100	100	100	100	n/a	n/a
0.025	54.5	48.5	44.7	51.5	49.8	4.2	8.4
0.05	22.0	15.0	15.0	17.1	17.3	3.3	19.3
0.25	8.1	6.9	6.4	0	5.3	3.6	68.1
IC ₅₀ (ng/mL)	0.025	0.023	0.022	0.026			

Table SI2. Assessment of waveguide reproducibility by comparing matrix calibration curves on 3 batches over 3 days.

Concentration	Batch 2	Batch 3	Batch 4	Mean	sd	CV (%)
(ng/mL)	Normalised Signal (<i>n</i> = 4)	Normalised Signal (<i>n</i> = 4)	Normalised Signal (<i>n</i> = 4)			
0	100	100	100	100	n/a	n/a
0.02	62.3	61.7	62.3	62.1	0.3	0.56
0.04	36.0	35.6	33.9	35.2	1.1	3.2
0.08	5.5	16.3	8.3	10.0	5.6	55.9
IC ₅₀ (ng/mL)	0.026	0.022	0.025			

- a. Paralytic shellfish toxins can enter the food web risking human health.
- b. Coastal monitoring of *Alexandrium spp.* can provide an early warning of risk.
- c. A highly innovative sample preparation method was developed.
- d. A sensitive rapid, portable assay was developed for paralytic shellfish toxins in seawater.
- e. The test was proven to accurately detect toxin presence in algae an unprecedented cell density.

Accepted manuscript

Fig. 1.

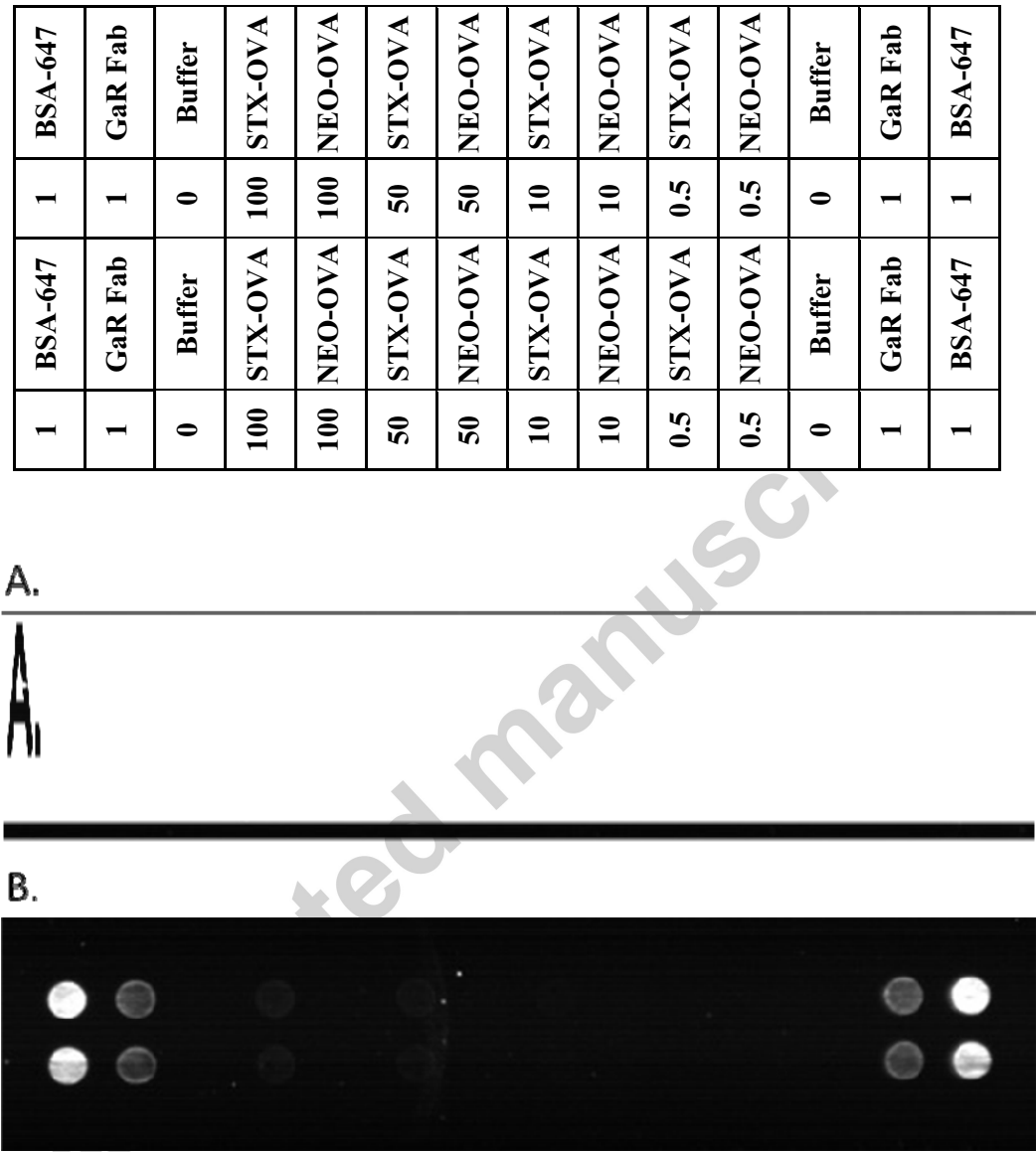


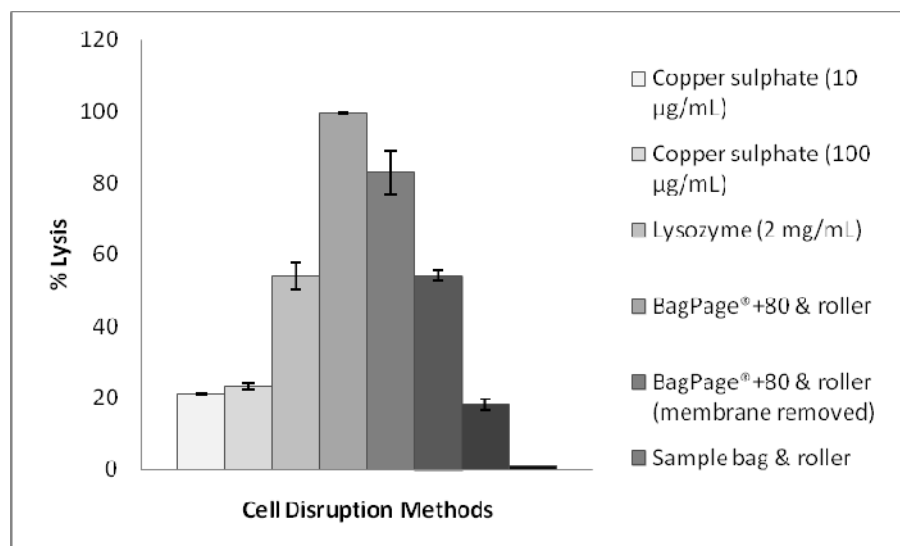
Fig. 2.

Fig. 3.

