



Method Article

Optimization of Loop-Mediated Isothermal Amplification (LAMP) reaction mixture for biosensor applications



Sümevra Vural Kaymaz, Meltem Elitas*

Faculty of Engineering and Natural Sciences, Sabanci University, Istanbul 34956, Turkey

ABSTRACT

Genetically Modified (GM) foods are becoming the future of agriculture on surviving global natural disasters and climate change by their enhanced production efficiency and improved functional properties. On the other hand, their adverse health and environmental effects, ample evidence on transgene leakage of Genetically Modified Organisms (GMOs) to crops have raised questions on their benefits and risks. Consequently, low-cost, reliable, rapid, and practical detection of GMOs have been important. GMO-detection platforms should be capable of stably storing detection reagents for long-delivery distances with varying ambient temperatures. In this study, we developed an event-specific, closed tube colorimetric GMO detection method based on Loop-Mediated Isothermal Amplification (LAMP) technique which can be integrated into GMO-detection platforms. The entire detection process optimized to 30 min and isothermally at 65 °C. The durability of the LAMP mixture in the test tubes showed that the LAMP reaction mixture, in which *Bst* polymerase and DNA sample was later included, yielded DNA amplicons for 3 days at room temperature, and for 6 days at 4 °C.

- Simple, stable, and cheap storage method of LAMP reaction mixture for GMO-detection technologies.
- GMO-detection platforms can stably store detection reagents for long-delivery distances with varying ambient temperatures.
- Any DNA sample can be used in the field or resource-limited setting by untrained personnel.

© 2021 The Author(s). Published by Elsevier B.V.

This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>)

ARTICLE INFO

Method name: Loop-Mediated Isothermal Amplification (LAMP) for genetically modified food detection

Keywords: Loop-Mediated Isothermal Amplification (LAMP), Genetically Modified Foods (GM foods), Resource-limited settings, Storage

Article history: Received 14 September 2020; Accepted 15 February 2021; Available online 19 February 2021

* Corresponding author.

E-mail address: melitas@sabanciuniv.edu (M. Elitas).

Specifications table

Subject area:	Biosensors
More specific subject area:	Storage of reaction reagents in biosensors
Method name:	Loop-Mediated Isothermal Amplification (LAMP) for genetically modified food detection
Name and reference of original report:	10.1016/j.bios.2019.111409
Resource availability:	Yes

Introduction

Genetically Modified Organisms (GMOs) products have been widespread and the increase in their production will undoubtedly continue with the discussion on the benefits and risks of Genetically Modified (GM) foods. Particularly, a variety of hypothesized potential adverse health and environmental effects of GM foods have made their detection methods essential. Along this line, GMO-detection platforms should be affordable, hand-held, easy-to-use, field-deployed, and their test results should be real-time [13]. Still, most of the detection methods for GM products relies on well-established laboratory techniques those require trained personnel, expensive laboratory settings and interpretation of obtained results to conclude the tests. Polymerase Chain Reaction (PCR), Real-Time PCR (RT-PCR), or digital PCR (dPCR) are among the broadly performed nucleic acid amplification-based detection tests in the laboratories [15,19,26]. Although they have high sensitivity, they demand expensive and specialized equipment that can be run by a specialist. PCR-based tests require longer test run time [15,19,26]. Moreover, their readout cannot be visualized by naked eye in real time.

In the past two decades, limitations of PCR have been overcome by isothermal nucleic acid amplification techniques [3,11,28]. These methods do not require cyclic temperature changes, instead it can be conducted at a constant temperature using simple instruments – e.g., a conventional water bath or heat block – which drastically reduces equipment costs [4,5,11]. Hence, fast and inexpensive detection methods can be used in comparison to PCR methods. Several isothermal exponential amplification methods have been developed, among them Loop-Mediated Isothermal Amplification (LAMP) exponentially amplifies nucleic acid template using 4–6 primers and DNA polymerase and provides 10⁹ copies of target DNA at 60–65°C in less than an hour [20]. LAMP provides high selectivity, specificity, efficiency and rapidity; therefore, it becomes an adequate DNA amplification method for GMO-detection platforms [4,18,21,22,24,30]. DNA-detection technologies mostly focused on providing rapid, sensitive, and accurate detection in cheap, reusable, portable devices those Rese [12]–[14].

However, when they include reaction mixtures, their storage and delivery, which directly contributes their performance and affordability, are also become an important technical challenge to solve. Since, these novel DNA-detection platforms use thermally unstable reagents for the enzymatic reactions to amplify nucleic acids, they must be stored and transported in refrigerated conditions [6,9,16,17,27]. This common limitation of DNA-detection platforms restricts their usage in resource limited settings, where reliable cold-chain equipment is not available [2,6,23]. Therefore, the reagents, those will be used in the DNA-detection platforms, should be carefully considered. Moreover, stability of the reagents inside the DNA-detection biosensors might contribute to determine their shelf life.

In this study, we exhibited the durability of the LAMP reaction mixture in the test tubes, excluded from *Bst* polymerase enzyme and DNA sample, for 3 days at room temperature (RT) and 6 days at 4°C. Our results might enable rapid detection of nucleic acids in resource-limited fields where transportation and storage conditions are not favorable. Herein, we showed the readout of the most convenient LAMP reaction mixture and its storage conditions using the lectin gene of soybean (*Glycine max*). Lectin gene has been widely used as a species control for all soybean samples. The P35S primer pair is used to detect transition site from the Cauliflower Mosaic Virus 35S promoter sequence in Roundup Ready Soybean (RRS, containing 10% gts40–3–2/RRS). We used RRS with P35S primer set to present specificity of the LAMP reactions.

Experimental section

DNA extraction and primer sequences

Soybean, MON89788 (100% GMO) is certified reference material (CRM), was obtained from the American Oil Chemists' Society (Boulder, Urbana, USA). Certified reference materials for Roundup Ready Soybean (RRS, gts40-3-2, 10% GMO) was obtained from Sigma Aldrich (St. Louis, MO, USA). Overall DNA extraction completed from 200 mg of soybean sample using the Macherey-Nagel NucleoSpin Plant Purification II kit according to the manufacturers' instructions. At the final step, 40 ml Elution Buffer was used 2 times to recover DNA for obtaining the highest concentration. DNA yield and purity were evaluated by UV Spectrophotometry at 230, 260 and 280 nm using a NanoDrop 2000c instrument (Thermo Scientific, Wilmington, DE, USA). DNA integrity was detected by agarose gel electrophoresis, in which 400–1200 ng/25 μ l DNA samples were separated on 1% agarose gels containing GelRed nucleic acid stain (Biotium, Hayward, CA, USA) in 0.5 Tris/Borate/EDTA (TBE) buffer. Primer sequences were reported in our previous work. The oligonucleotides were obtained from Oligomer (Ankara, TURKEY) [10,13].

LAMP reaction

Each LAMP reaction was performed in a final volume of 25 μ l containing 8 U *Bst* DNA polymerase (large fragment; New England Biolabs) in 1X ThermoPol Reaction buffer (20 mM Tris-HCl, 10 mM KCl, 10 mM $(\text{NH}_4)_2\text{SO}_4$, 2 mM MgSO_4 and 0.1% Triton X-100; New England Biolabs) supplemented with 6 mM MgSO_4 , 1 M betaine and 1.4 mM of each dNTP, 2 μ l DNA template, 6.5 μ l double-distilled water (ddH_2O), and 1.6 μ l, 1X primer mix. Each reaction was incubated at 65 °C for 30 min in a heat block. LAMP results were determined according to color change provided by Hydroxy Naphthol Blue (HNB), which was added in 120 μ M to each test reaction before the incubation.

Shelf life of LAMP reaction reagents

The preparation of the reactive mixtures was kept both at room temperature (25°C) and in the refrigerator at 4 °C for 7 days to test their shelf life. LAMP reagents were mixed in a PCR tube and stored wrapped in aluminum foil for preventing light transmittance. Next, efficiency of the LAMP reaction was tested by introducing the sample DNA and *Bst* polymerase enzyme at 65 °C for 30 min.

Detection and quantification of the LAMP product

When the reaction incubation completed, the color of each reaction was assessed by naked eye where negative results were differentiated from positive results by a color change from violet to sky blue. LAMP reaction products were also resolved electrophoretically (1% agarose in 0.5X TBE buffer) and visualized using BioRad nucleic acid stain in BioRad Gel Doc EZ System (Berkeley, CA, USA).

Data quantification was performed using ImageJ software (Version 2.0 National Institutes of Health, Rockville, MD, USA). Intensity of the bands in the images was acquired from the agarose gel images by removing the backgrounds of the DNA bands and defining a rectangular region of interest. Then, lane profile plots were drawn based on the measured areas (Analyze-Gels-Plot Lanes). The obtained data was analyzed and presented using GraphPad Prism software (Version 5). Student's *t*-test was used to determine statistical significance of changes in band intensities. Figures show the data as mean $\pm$ standard deviation. The experiments were performed in triplicate.

Results and discussion

Currently, widespread application of molecular biology techniques in research, diagnosis, and education is costly and often unaffordable for many laboratories and schools in resource-poor countries [2]. Particularly, in case of colorimetric tests, it is important to know the stability of the indicator to facilitate the use of detection platforms in the field [1]. In this study, we reported a LAMP

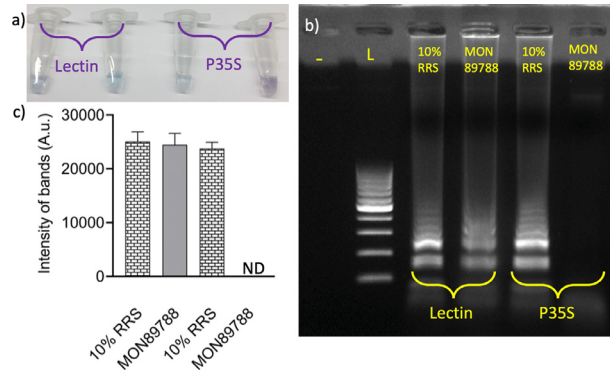


Fig. 2. Specificity of LAMP-based detection. LAMP reactions using lectin and P35S primer sets with RSS (10% gts40–3–2/RSS soybean) and MON89788 (100% GMO) template DNA. **(a)** Colorimetric results. **(b)** Gel electrophoresis results with the order of the samples shown in a. 1-kbp DNA ladder is used. **(c)** Intensity measurements of the amplified DNA fragments following the order of the samples shown in b (For interpretation of the references to color in this figure, the reader is referred to the web version of this article.). Intensity measurements of the DNA bands show mean $\pm$ standard deviation of three independent experiments. ND: Not detected.

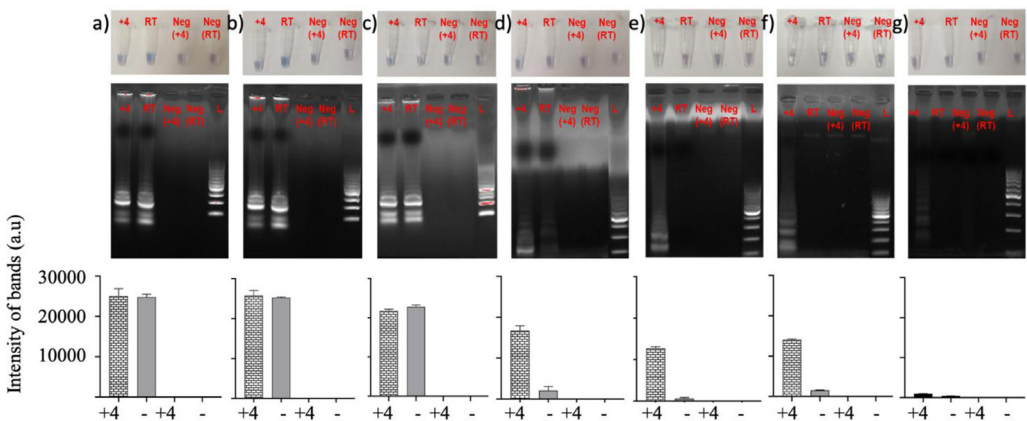


Fig. 3. Testing stability of the LAMP reaction mixtures (RM) at RT (25°C) and at 4°C for 7 days. LAMP reactions using the lectin primer set with 100 ng of MON89788 template DNA were incubated for **(a)** 1, **(b)** 2, **(c)** 3, **(d)** 4, **(e)** 5, **(f)** 6, **(g)** 7 days either at RT or 4°C fridge. Neg: Negative control, LAMP reaction without DNA template. L: 1-kbp DNA ladder. Results were presented by gel electrophoresis (1% agarose gel), colorimetric images (images of the reaction tubes) and the intensity measurements of the DNA bands.

storage conditions, their repeated freezing and thawing can lead to loss of reactivity in their long-term storage. For a regular LAMP assay, DNA polymerase, primers, and dNTPs need to be stored at -20°C . Herein, we verified that LAMP reaction mixture without DNA template and polymerase enzyme, can be stored at 4°C for 6 days and at RT for 3 days. The results were verified by DNA amplicons on agarose gel electrophoresis, Fig. 3.

In this study, we showed the readout of the most convenient LAMP reaction mixture and its storage conditions using the lectin gene of soybean. Lectin gene has been widely used as a housekeeping gene to confirm soybean species. The storage of the reaction mixture is not target gene or species specific. When different genes from different species are identified using LAMP protocol, storage of their reaction mixture will remain the same. Our results suggest that LAMP reaction mixture (RM 5) without DNA template and polymerase enzyme can be delivered in a biosensor or a PCR test tube at 4°C for 6 days and at RT for 3 days. The end-user can only add the enzyme and

Table 1
Lectin and P35S primer sets used in this study [10].

Lectin-F3:	ACAACTTTTGAAAAGTACCCAAT
Lectin-B3:	GAAC TTGTTCCAGCTGAAAG
Lectin-FIP:	CATTGCTTTGCTTCAGCTAAATTGCTGCTAGTATAAATAGGGGCATG
Lectin-BIP:	AACCCAGAATGTGGTTGTATCTCTCAGTTTCCGCTGAGTTTGC
P35S_LampF:	GTCTTCAAAGCAAGTGTTTGGATAGTGGGATTGTGCG
P35S_LampR:	TTCCACGATGCTCCTCGTTTCTCTGCCGACAGTGG
P35S_LoopF:	TCCACTGACGTAAGGG
P35S_LoopR:	GGGGTCCATCTTTGGG
P35S_DisplF:	AGGAAGGGTCTTGCG
P35S_DisplR:	ATAAAGGAAAGGCCATCG

Table 2
The tested concentrations of reagents for optimization of LAMP readout.

Reagents	RM1	RM2	RM3	RM4	RM5
1X Isothermal Amplification Buffer (μl)	2.5	2.5	2.5	2.5	2.5
MgSO ₄ (mM)	4	8	6	8	6
10x Primers (μl)	2.5	2.5	2.5	2.5	2.5
10 mM dNTPs (μl)	3.5	3.5	3.5	3.5	3.5
Betaine (μl)	5	5	5	5	5
HNB (μM)	120	120	240	240	120
<i>Bst</i> 2.0 (μl)	1	1	1	1	1
Template DNA (ML)	2	2	2	2	2
Double Distilled Water (μl)	7.5	6.5	6.5	6.5	7.5
Total volume (μl)	25	25	25	25	25

sample DNA into the reaction which can be practical for an untrained person. Besides, LAMP assays require 60–65°C temperature by a heating block or water bath, which is mostly available in many resource-limited settings and does not require expertise for their operation. Therefore, our results might provide a solution for LAMP-based colorimetric DNA-detection biosensors containing premixed reaction reagents.

Conclusion

Stability of reaction mixtures is an important parameter for the newly designed DNA detection platforms which includes premixed reaction reagents. In this study, a LAMP based colorimetric GMO detection method was optimized for testing the durability of LAMP reaction mixtures for 7 days both at RT and 4°C. First, LAMP reaction mixture was prepared without DNA template and *Bst* polymerase and stored either at RT or 4°C. Next, DNA template and *Bst* polymerase was added into the reaction mixture and the LAMP reaction was performed at 65 °C for 30 min. Afterwards, color change was visualized, the blue color indicated the existence of GMO, positive reactions. Finally, verification of the results was performed by gel electrophoresis. Our method presents that LAMP reaction can be performed when the reaction mixture without polymerase enzyme and DNA template is stored for 3 days at RT and 6 days at 4°C.

Table 1

Declaration of Competing Interest

The authors declare no competing interests.

Acknowledgments

This work was supported by Scientific and Technological Research Council of Turkey (TUBITAK) under grant 118Z217.

References

- [1] B.S. Antharavally, K.A. Mallia, P. Rangaraj, P. Haney, P.A. Bell, Quantitation of proteins using a dye-metal-based colorimetric protein assay, *Anal. Biochem.* 385 (2) (2009) 342–345, doi:[10.1016/j.ab.2008.11.024](https://doi.org/10.1016/j.ab.2008.11.024).
- [2] S. Bhadra, A. Pothukuchy, R. Shroff, A.W. Cole, M. Byrom, J.W. Ellefson, J.D. Gollihar, A.D. Ellington, Cellular reagents for diagnostics and synthetic biology, *PLOS ONE* 13 (8) (2018) 1–24, doi:[10.1371/journal.pone.0201681](https://doi.org/10.1371/journal.pone.0201681).
- [3] O.L. Bodulev, I.Y. Sakharov, Isothermal nucleic acid amplification techniques and their use in bioanalysis, *Biochemistry (Mosc.)* 85 (2) (2020) 147–166, doi:[10.1134/S0006297920020030](https://doi.org/10.1134/S0006297920020030).
- [4] T. Deng, W. Huang, R. Xing, J. Zhang, J. Ren, Y. Ge, Y. Chen, Establishment and application of a loop-mediated isothermal amplification method with double-stranded displacement probes to quantify the genetically modified rice M12 event, *Eur. Food Res. Technol.* 246 (3) (2020) 631–641, doi:[10.1007/s00217-020-03430-1](https://doi.org/10.1007/s00217-020-03430-1).
- [5] F. Carballo, B. Leticia, (2017). Low-cost point-of-care diagnostic devices for low resource settings. TDX (Tesis Doctorals En Xarxa). <http://www.tdx.cat/handle/10803/401780>
- [6] A. Germini, S. Rossi, A. Zanetti, R. Corradini, C. Fogher, R. Marchelli, Development of a peptide nucleic acid array platform for the detection of genetically modified organisms in food, *J. Agric. Food Chem.* 53 (10) (2005) 3958–3962, doi:[10.1021/jf050016e](https://doi.org/10.1021/jf050016e).
- [7] P. Gill, A. Ghaemi, Nucleic acid isothermal amplification technologies - a review, *Nucleosides, Nucleotides Nucleic Acids* 27 (3) (2008) 224–243, doi:[10.1080/15257770701845204](https://doi.org/10.1080/15257770701845204).
- [8] M. Goto, E. Honda, A. Ogura, A. Nomoto, K.I. Hanaki, Colorimetric detection of loop-mediated isothermal amplification reaction by using hydroxy naphthol blue, *BioTechniques* 46 (3) (2009) 167–172, doi:[10.2144/000113072](https://doi.org/10.2144/000113072).
- [9] L. Gu, W. Yan, L. Liu, S. Wang, X. Zhang, M. Lyu, Research progress on rolling circle amplification (Rca)-based biomedical sensing, *Pharmaceuticals* 11 (2) (2018) 1–19, doi:[10.3390/ph11020035](https://doi.org/10.3390/ph11020035).
- [10] X. Guan, J. Guo, P. Shen, L. Yang, D. Zhang, Visual and rapid detection of two genetically modified soybean events using loop-mediated isothermal amplification method, *Food Anal. Methods* 3 (4) (2010) 313–320, doi:[10.1007/s12161-010-9132-x](https://doi.org/10.1007/s12161-010-9132-x).
- [11] L. Guo, B. Lu, Q. Dong, Y. Tang, Y. Du, B. Li, O. Mukama, J. de Dieu Habimana, X. Meng, Y. Ting, F. Songwe, A. Al Farga, S. Mugisha, P. Rwibasira, Z. Li, L. Zeng, One-tube smart genetic testing via coupling isothermal amplification and three-way nucleic acid circuit to glucometers, *Anal. Biochem.* 1106 (December 2019) (2020) 113762, doi:[10.1016/j.ab.2020.113762](https://doi.org/10.1016/j.ab.2020.113762).
- [12] P. Herman, I.V.J. Taverniers, M. Fraiture, P. Herman, I. Taverniers, M.De Loose, D. Deforce, N.H. Roosens, Current and New Approaches in GMO Detection : Challenges and Solutions, *BioMed Res. Int.* 2015 (2015) (Table 1), doi:[10.1155/2015/392872](https://doi.org/10.1155/2015/392872).
- [13] D. Kaygusuz, S. Vural, A.Ö. Aytekin, S.J. Lucas, M. Elitas, DaimonDNA: a portable, low-cost loop-mediated isothermal amplification platform for naked-eye detection of genetically modified organisms in resource-limited settings, *Biosens. Bioelectron.* 141 (March) (2019) 1–6, doi:[10.1016/j.bios.2019.111409](https://doi.org/10.1016/j.bios.2019.111409).
- [14] M.S.P. López, C.L. Manzanera-Palenzuela, B. López-Ruiz, Biosensors for GMO testing: nearly 25 years of research, *Crit. Rev. Anal. Chem.* 48 (5) (2018) 391–405, doi:[10.1080/10408347.2018.1442708](https://doi.org/10.1080/10408347.2018.1442708).
- [15] X. Liu, J. Jiang, K. Zhang, E. Long, J. Cui, M. Zhu, Y. An, J. Zhang, Z. Liu, Z. Lin, X. Li, J. Chen, Q. Cao, J. Li, X. Wu, D. Wang, H. Lin, Localization and diagnosis framework for pediatric cataracts based on slit-lamp images using deep features of a convolutional neural network, *PLOS ONE* 12 (3) (2017) 1–18, doi:[10.1371/journal.pone.0168606](https://doi.org/10.1371/journal.pone.0168606).
- [16] S. Lutz, P. Weber, M. Focke, B. Faltin, J. Hoffmann, C. Müller, D. Mark, G. Roth, P. Munday, N. Armes, O. Piepenburg, R. Zengerle, F. Von Stetten, Microfluidic lab-on-a-foil for nucleic acid analysis based on isothermal recombinase polymerase amplification (RPA), *Lab Chip* 10 (7) (2010) 887–893, doi:[10.1039/b921140c](https://doi.org/10.1039/b921140c).
- [17] J. Mano, N. Shigemitsu, S. Futo, H. Akiyama, R. Teshima, A. Hino, S. Furui, K. Kitta, Real-time PCR array as a universal platform for the detection of genetically modified crops and its application in identifying unapproved genetically modified crops in Japan, *J. Agric. Food Chem.* 57 (1) (2009) 26–37, doi:[10.1021/jf802551h](https://doi.org/10.1021/jf802551h).
- [18] Y. Mori, T. Notomi, Loop-mediated isothermal amplification (LAMP): expansion of its practical application as a tool to achieve universal health coverage, *J. Infect. Chemother.* 26 (1) (2020) 13–17, doi:[10.1016/j.jiac.2019.07.020](https://doi.org/10.1016/j.jiac.2019.07.020).
- [19] I. Nazarenko, B. Lowe, M. Darfler, P. Ikonomi, D. Schuster, A. Rashtchian, Multiplex quantitative PCR using self-quenched primers labeled with a single fluorophore, *Nucleic Acids Res.* 30 (9) (2002) e37, doi:[10.1093/nar/30.9.e37](https://doi.org/10.1093/nar/30.9.e37).
- [20] T. Notomi, H. Okayama, H. Masubuchi, W. Yonekawa, A. Keiko, H.T. Nobuyuki, Loop-mediated isothermal amplification of DNA, *Nucleic Acids Res.* 28 (2000) e63 doi:[10.1093/nar/28.12.e63](https://doi.org/10.1093/nar/28.12.e63).
- [21] S. Panno, S. Matić, A. Tiberini, A.G. Caruso, P. Bella, L. Torta, R. Stassi, S. Davino, Loop mediated isothermal amplification: principles and applications in plant virology, *Plants* 9 (4) (2020) 1–28, doi:[10.3390/plants9040461](https://doi.org/10.3390/plants9040461).
- [22] J.C. Patel, N.W. Lucchi, P. Srivastava, J.T. Lin, R. Sug-Aram, S. Aruncharus, P.K. Bharti, M.M. Shukla, K. Congpuong, W. Satimai, N. Singh, V. Udhayakumar, S.R. Meshnick, Field evaluation of a real-time fluorescence loop-mediated isothermal amplification assay, realamp, for the diagnosis of Malaria in Thailand and India, *J. Infect. Dis.* 210 (8) (2014) 1180–1187, doi:[10.1093/infdis/jiu252](https://doi.org/10.1093/infdis/jiu252).
- [23] B. Seise, S. Pollok, C. Seyboldt, K. Weber, Dry-reagent-based PCR as a novel tool for the rapid detection of *Clostridium* spp, *J. Med. Microbiol.* 62 (PART10) (2013) 1588–1591, doi:[10.1099/jmm.0.060061-0](https://doi.org/10.1099/jmm.0.060061-0).
- [24] M. Singh, D. Pal, P. Sood, G. Randhawa, Construct-specific loop-mediated isothermal amplification: rapid detection of genetically modified crops with insect resistance or herbicide tolerance, *J. AOAC Int.* (2020) March, 1–10, doi:[10.1093/jaoacint/qsaa043](https://doi.org/10.1093/jaoacint/qsaa043).
- [25] L.D. Thiessen, T.M. Neill, W.F. Mahaffee, Development of a quantitative loop-mediated isothermal amplification assay for the field detection of erysiphe necator, in: *PeerJ*, 6, 2018, p. e4639, doi:[10.7717/peerj.4639](https://doi.org/10.7717/peerj.4639).
- [26] W. Tian, P. Li, W. He, C. Liu, Z. Li, Rolling circle extension-actuated loop-mediated isothermal amplification (RCA-LAMP) for ultrasensitive detection of microRNAs, *Biosens. Bioelectron.* 128 (December 2018) (2019) 17–22, doi:[10.1016/j.bios.2018.12.041](https://doi.org/10.1016/j.bios.2018.12.041).
- [27] H. Wang, Z. Ma, J. Qin, Z. Shen, Q. Liu, X. Chen, H. Wang, Z. An, W. Liu, M. Li, A versatile loop-mediated isothermal amplification microchip platform for streptococcus pneumoniae and mycoplasma pneumoniae testing at the point of care, *Biosens. Bioelectron.* 126 (2019) 373–380, doi:[10.1016/j.bios.2018.11.011](https://doi.org/10.1016/j.bios.2018.11.011).

- [28] P.L. White, A. Alanio, M. Cruciani, R. Gorton, L. Millon, V. Rickerts, R.A. Barnes, J.P. Donnelly, J. Loeffler, Nucleic acid tools for invasive fungal disease diagnosis, *Curr. Fungal Infect. Rep.* 14 (1) (2020) 76–88, doi:[10.1007/s12281-020-00374-7](https://doi.org/10.1007/s12281-020-00374-7).
- [29] M. Yan, W. Li, Z. Zhou, H. Peng, Z. Luo, L. Xu, Direct detection of various pathogens by loop-mediated isothermal amplification assays on bacterial culture and bacterial colony, *Microb. Pathog.* 102 (2017) 1–7, doi:[10.1016/j.micpath.2016.10.025](https://doi.org/10.1016/j.micpath.2016.10.025).
- [30] Z. Yun, L. Peng, S. Huo, X. Wu, J. Zhou, H. Zhang, W. Li, C. Qi, G. Huang, Auto-microfluidic thin-film chip for genetically modified maize detection, *Food Control* 80 (2017) 360–365, doi:[10.1016/j.foodcont.2017.03.035](https://doi.org/10.1016/j.foodcont.2017.03.035).
- [31] M. Zhang, J. Ye, J.song He, F. Zhang, J. Ping, C. Qian, J. Wu, J.J. Li, C. Xiong, Y. Liu, J.S. Liang, X.W. Zhou, Visual detection for nucleic acid-based techniques as potential on-site detection methods. A review, *Anal. Chim. Acta* 7 (DECEMBER2016) (2020) 1–15, doi:[10.3389/fpls.2016.01956](https://doi.org/10.3389/fpls.2016.01956).